





EUROLINE Infectious Serology Profiles

Products available

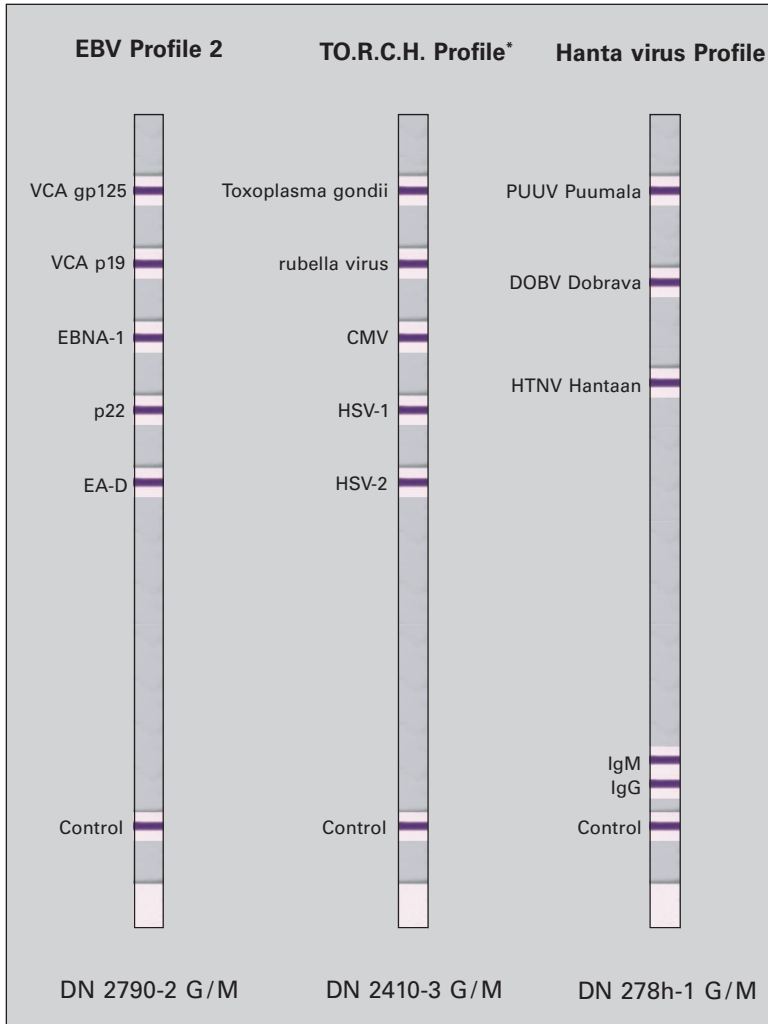

*Currently not available as IVD in the European Union.

Examples of incubated strips

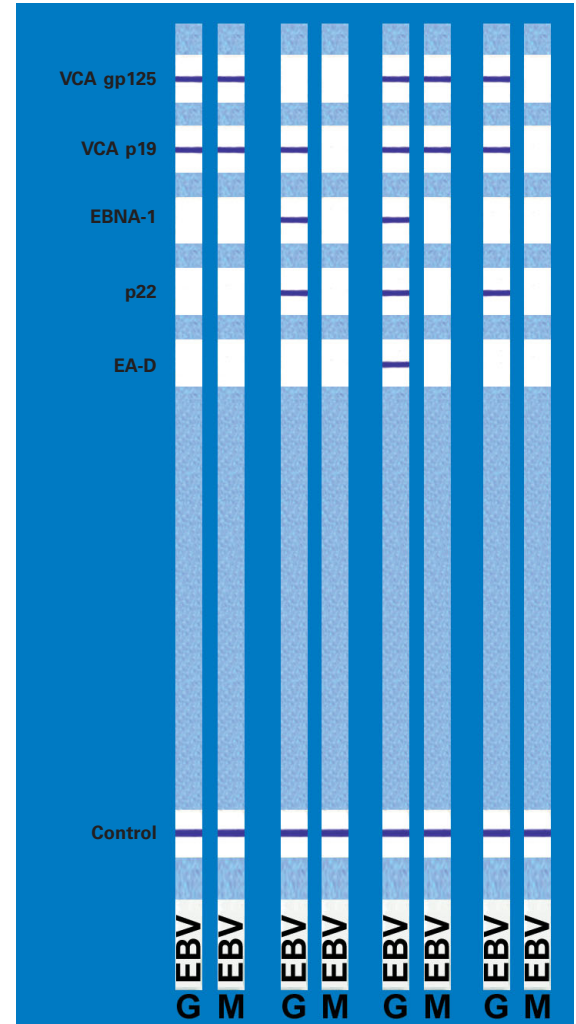

Technical data

Applications	Determination of antibodies against infectious agent antigens in the serological diagnosis of infectious diseases.
Antigen substrates	Membrane chips printed with individual lines of purified, biochemically characterized antigens. Each antigen is coated onto a separate membrane chip, providing optimal efficiency of antibody detection for each antigen.
Sample dilution	Serum or plasma; 1:101 in buffer.
Test procedure	60 min/60 min/20 min (DN 2790); 30 min/30 min/10 min (DN 2410); 30 min/30 min/10 min (DN 278h-1 G/M), room temperature.
Automation options	Compatible with all commercial blot processing systems. Computer-based evaluation of results with EUROLineScan program.
Kit format	16 membrane strips. Kits include all necessary reagents.